



WalkMed® 350VL™ Ambulatory Infusion Pump

High Precision, High Reliability

The delivery system in WalkMed ambulatory infusion pumps is highly regarded for its trouble-free performance and its ability to deliver accurately, even at low infusion rates.

The WalkMed 350VL is ideal for chemotherapy. The new features of the WalkMed 350VL build on the reputation established by the widely utilized WalkMed 350 chemotherapy pump.

Operational Highlights

- ◆ Easy 4-step Programming
- ◆ 0.10 – 30.0 mL/hr Delivery Rates
- ◆ Compact & Durable

Safety & Convenience

- ◆ Keypad Lock-out for Patient Protection
- ◆ Auto-Occlusion Alarm Release for Sensing & Silencing Alarms
- ◆ Programmable Total Volume Limit & Delivery Tracking

Need help with a WalkMed brand product? Call our 24-Hour Support Line: (866) 484-5601.

CAUTION: Federal (USA) law restricts this device to sale by, or on the order of a physician.




The WalkMed® 350VL™

provides clinicians with the dependability they have come to expect from a company with over 20 years of ambulatory infusion experience and provides the features that a growing and dynamic market demand.

The WalkMed brand is recognized by the medical community as simple and easy to use. The new WalkMed 350VL has been designed to improve user friendliness and enhance the therapy experience for patients.

Easy to Use Keypad

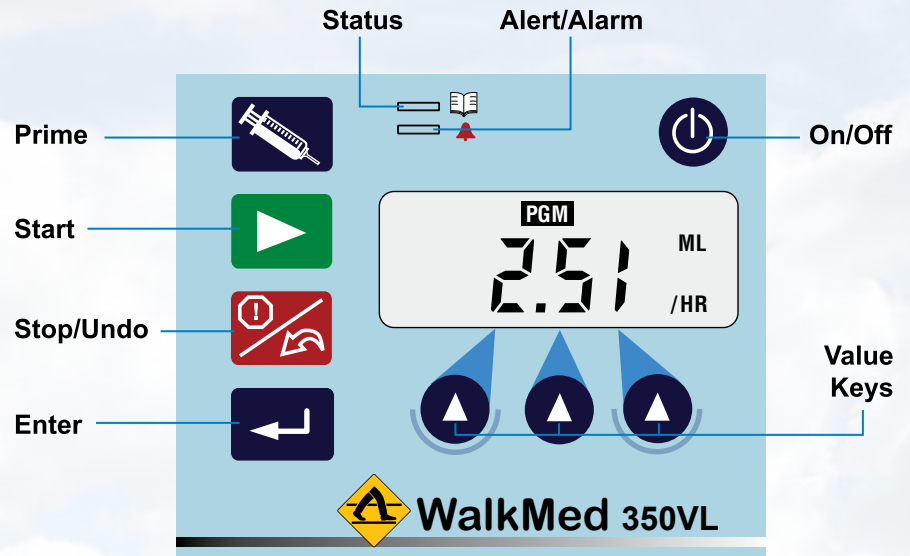

Actual Size

Specifications

Dimensions	4.4 x 4.0 x 1.5 in (includes front cover)
Weight	12 oz (includes battery and front cover)
Accuracy	+/- 5%
Infusion Mode	Continuous
Flow Rate	0.10 – 19.99 mL/hr (0.01 mL resolution) 20.0 – 30.0 mL/hr (0.10 mL resolution)
Total Volume Limit	1 – 1,999 mL
Reservoir Volume	65, 150, 250 mL bags (spike sets available)
Power	9 volt battery (standard)
Battery Capacity	450 mL @ 1 mL/hr - 18.7 Days 650 mL @ 10 mL/hr - 2.7 Days 500 mL @ 30 mL/hr - 0.7 Days
Alarms/Alerts	System Idle, Downstream Occlusion, Low Battery Depleted Battery, Under-Delivery, Over-Delivery, Not Delivering, Data Entry Error, Clamp Bar Open, System Malfunction, End of Infusion Notice, KVO/TVL Reached

New Features

- Patient Lock-Out
- Easier Menu System
- Increased Infusion Rate
- Total Volume Limit
- Volume Delivered Tracking
- Auto-Occlusion Alarm Release
- 4 Hour KVO Protocol
- End of Infusion Notice
- Symbolic Keypad
- Battery Retention Cover
- Snap On/Off Reservoir Covers

